

Lab #9 Biochemical activities of bacteria

To perform some biochemical tests which are routinely used in the identification of the different types of bacteria

Tests to perform

- H₂S production in SIM
- Indole
- Motility
- Lactose ,Glucose fermentation
- Starch hydrolysis
- Lipids

Biochemical Tests:


- They are done to measure the presence or absence of specific enzyme.
- It help identifying bacterial strains.
- Type of biochemical reaction act as a (thumbprint) bas on :
 - 1- each species has a different DNA
 - 2- different protein enzymes
 - 3- different and unique sets of biochemical reaction

Lipid hydrolysis:

- Some bacteria produce enzyme that hydrolyzes triglycerides into free fatty acids and glycerol which can then be taken up by the cell and metabolized more in glycolysis.
- **lipase** could be detected by growing the bacteria on **neutral red fat agar** - Neutral red (pH indicator that turns bright red at a low pH)

Lipid hydrolysis:

Positive
for lipid
hydrolysis

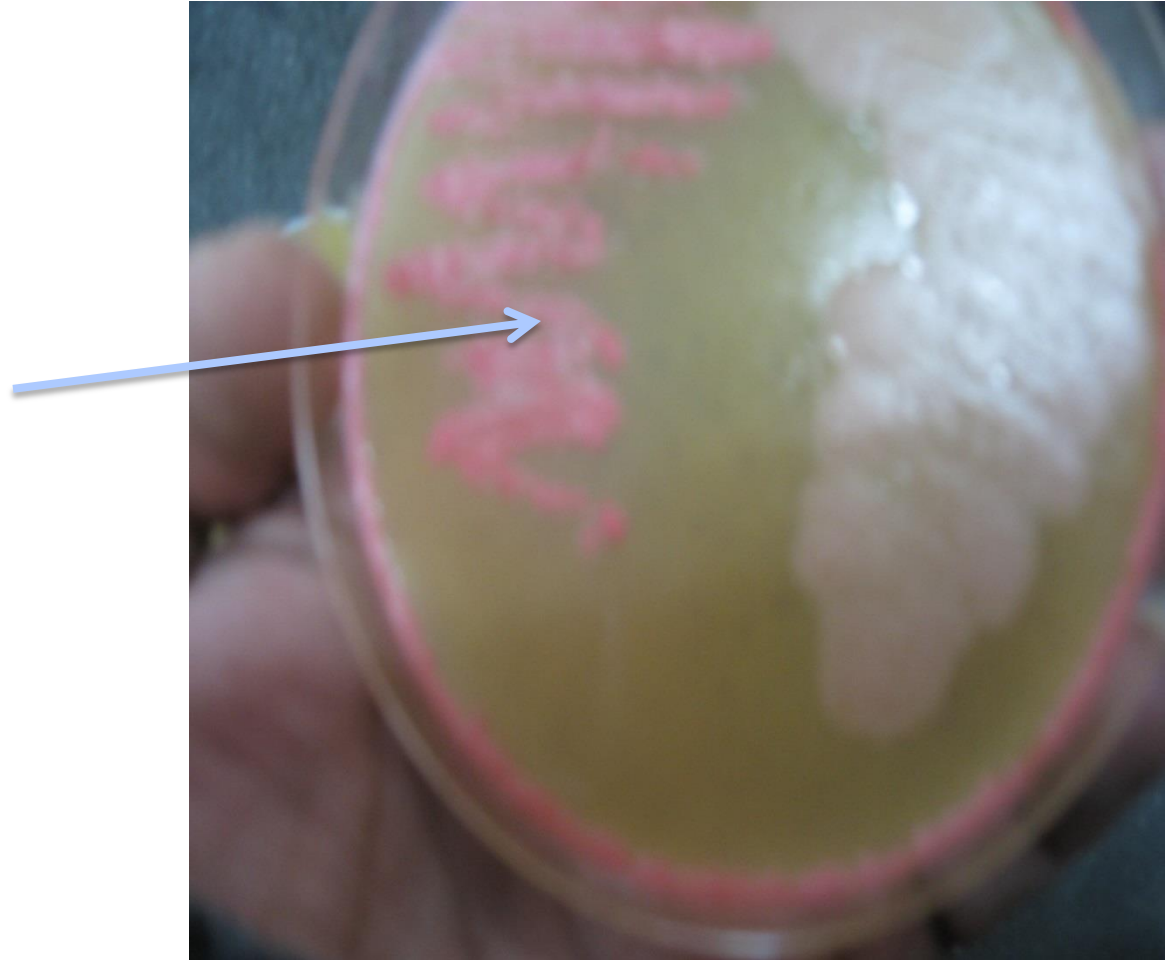

Gelatin liquefaction:

- Some bacteria can use gelatin as a source of carbon.
- **Protease** break the peptide bond between amino acids.
- Gelatin is semi-solid at room temperature. and melt at higher temperature such as 37C.
- In order to take the results, the tubes must be placed in the refrigerator for 5-10 minutes.

Gelatin Liquefaction

Can bacteria use gelatin as source of carbon ?

Negative result:-

Solid in ice bath or refrigerator after 5 minutes.



Positive result:

Liquid in ice bath or refrigerator after 5 minutes.

Stab straight to bottom with an inoculating needle. Incubate at 37C for 24 hr . Put in ice bath or refrigerator for 5 minutes, read results.

TSI

Purpose

Triple sugar iron (TSI) agar is a differential medium used in determining carbohydrate fermentation, H_2S production and Gas from carbohydrate metabolism.

- Bacteria can metabolize carbohydrates aerobically (with oxygen) or fermentatively (without oxygen). TSI differentiates bacteria based on their fermentation of lactose, glucose and sucrose and on the production of hydrogen sulfide.

- **Sugar fermentations**
- **Gas production**
- **Hydrogen sulfide production**

Pancreatic digest of casein	10.0 g
Peptic digest of animal tissue	10.0 g
Glucose	1.0 g
Lactose	10.0 g
Sucrose	10.0 g
Ferrous sulfate or ferrous ammonium sulfate	0.2 g
NaCl	5.0 g
Sodium thiosulfate	0.3 g
Phenol red	0.024 g
Agar	13.0 g
Distilled water	1,000 mL

TSI: Triple Sugar Iron Test.

- Triple sugar:
 - 1% lactose, 1% sucrose and 0.1% glucose.
 - The indicator is phenol red
 - A/A or yellow/ yellow → the bacteria can ferment the three sugars.
 - K/K or (red/red)→ non fermenters bacteria.
- Sulfur reduction:
 - The indicator is iron.
- Gas production:
 - When Peptone is utilized aerobically it produce ammonia gas.


SIM: Sulfur, Indole and Motility test

- Sulfur reaction:
 - ❖ Ferrous ammonium sulfate and sodium thiosulfate (together serve as indicator for the production **(H₂S)**).
 - ❖ H₂S react with iron and form **ferric sulfide. (black precipitate)**.
- Indole test:
 - ❖ Some bacteria can hydrolyze Tryptophan and produce **Indole** by **Tryptophanase**.
 - ❖ The indicator is Kovac`s reagent (positive results → red color).
- Motility test (semi-solid media)


SIM

Procedure

1. Using a wire needle , inoculate test organism two-thirds into the medium with stab motion.
2. Examine tubes after incubation for motility and H₂S production.
3. Add 2-3 drops of Kovac's Reagent to each tube. Record as **indole positive** if a pink or red color appear, or as **indole negative** if there is no color change.

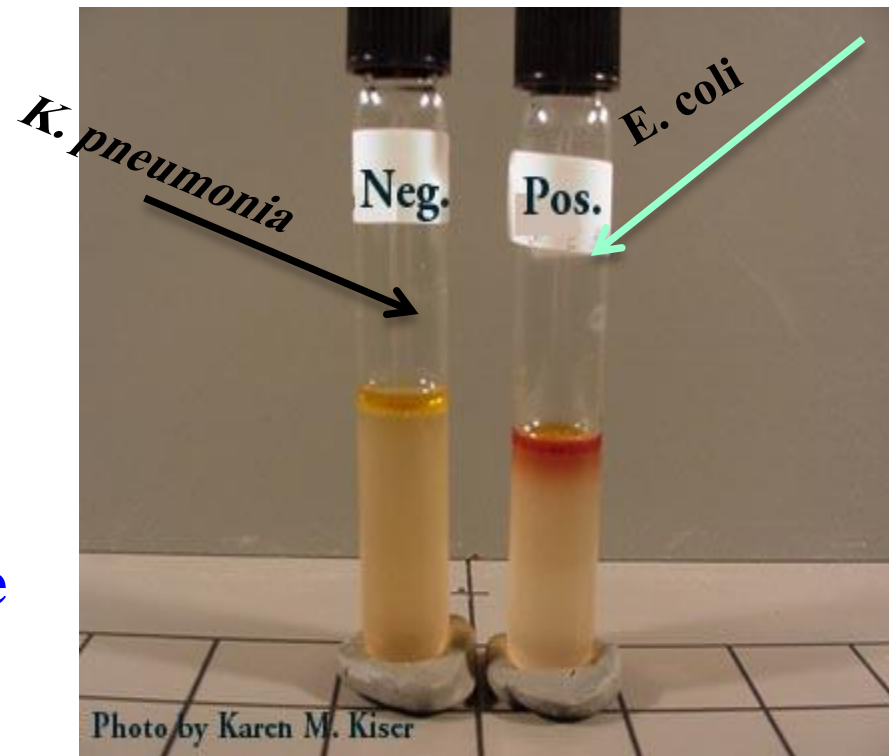


SIM: Sulfur, Indole and Motility test

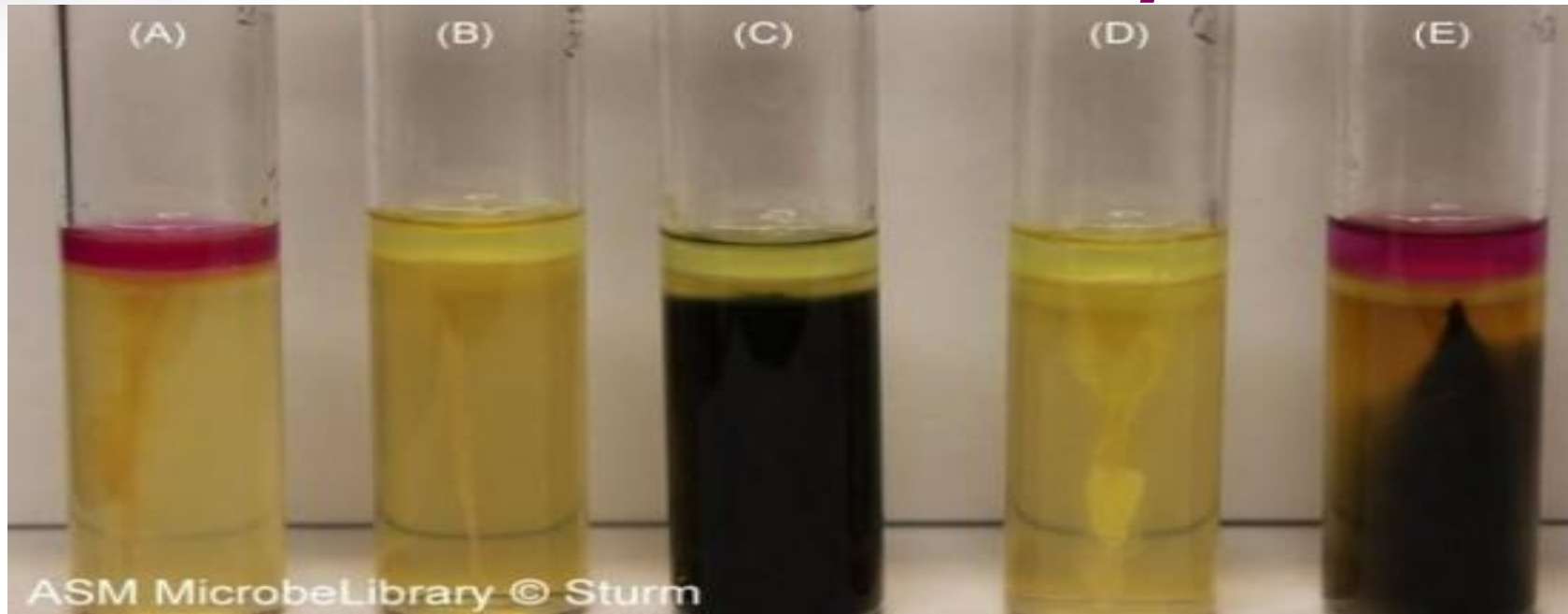

- **A: *Escherichia coli*** - Negative for H₂S, Positive for Indole, motile
- **B: *Klebsiella pneumoniae*** - Negative for H₂S, Negative for Indole, none motile
- **C: *Salmonella arizonae*** - Positive for H₂S, Negative for Indole, Positive for motility
- **D: *Enterobacter aerogenes*** - Negative for H₂S, Negative for Indole, Positive for motility
- **E: *Proteus vulgaris*** - Positive for H₂S, Positive for Indole, Positive for motility

Procedure

- 1- Use a straight inoculating loop to pickup an isolated colony.**
- 2- Inoculate the TSI slant by first stabbing the butt down to the bottom, withdraw the needle, and then streak the surface of the slant.**
- 3- Read results after incubation at 37°C for 18 to 24 h.**

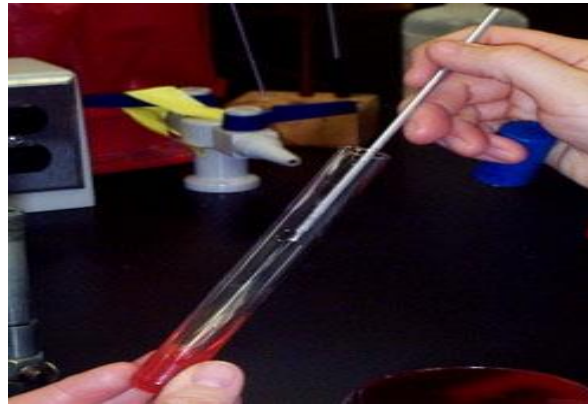


TSI: Triple Sugar Iron Test.



Starch Hydrolysis:

- Test the ability of an organism to produce certain coenzymes such as
- Starch addition make the media nutritive.
- NO color change occur in the medium when starch is hydrolyzed.
- IODINE (IKI) should be added to the plate AFTER incubation.
- If the color change to **black or blue** → no hydrolysis.
- If there was no change in medium color or a clear zone appeared → hydrolysis occurred

Starch Hydrolysis:

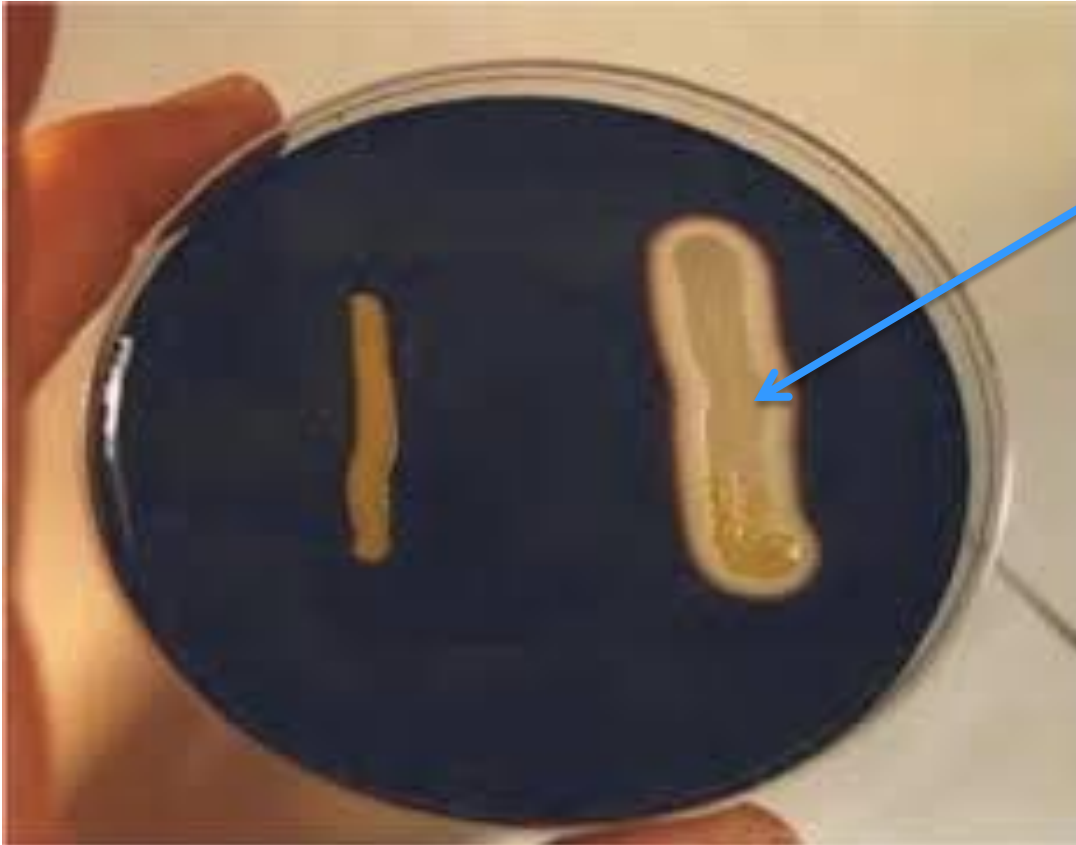

Positive for starch hydrolysis

Biochemical Tests:

Test	Test for	Indicator
SIM	Sulfur reduction / tryptophanase enzyme/ motility.	Iron for sulfur reduction / Kovac`s reagent for Indole.
TSI	Fermentation of three sugars / sulfur reduction/ gas production	Iron for sulfur reduction / phenol red → yellow (acid).
Starch hydrolysis	Starch hydrolysis by Amylase	Iodine
Lipid hydrolysis	Lipids hydrolysis by Lipases	Neutral red → red (acid)
Gelatin liquefaction	The presence of Proteases	-----
Citrate utilization	If the bacteria can use citrate as a source of carbon.	Bromothymo blue → blue (basic).
Urea utilization	The presence of Urease.	Phenol red→ Pink (basic)
Oxidation test	The presence of Cytochrome c oxidase.	Tetramethyl-p-phenyliamine → purple
Nitrate reduction	The presence of Nitroreductase	Alpha-naphtyl amine and sulfanic acid
Catalase test	Catalase enzyme	-----
Coagulase test	Coagulase enzyme	-----

Lab #10 Biochemical activities of bacteria (part 2)

Objectives

To perform some biochemical tests which are routinely used in the identification of the different types of bacteria

Tests to know

- 1- Citrate
- 2- Urea utilization
- 3- Oxidase
- 4- Nitrate reduction
- 5- Catalase
- 6- Coagulase
- 7-PAD

SIM

Citrate is an organic molecule that can be utilized by bacteria that produce the enzyme **citrase**. Citrase is produced by some bacteria such as *K.pneumoniae*, *E. aerogenes* but not by others like *E. coli*

Media and Reagent: Simmon's Citrate Agar. It has citrate as the only carbon source and pH indicator **bromothymol blue**


Method: Inoculate the slant and incubate at 37° C for 24-48 hours.

Expected results:

Positive test: Growth and color changes to blue

Negative test: No growth and color remains green


Citrate Test results

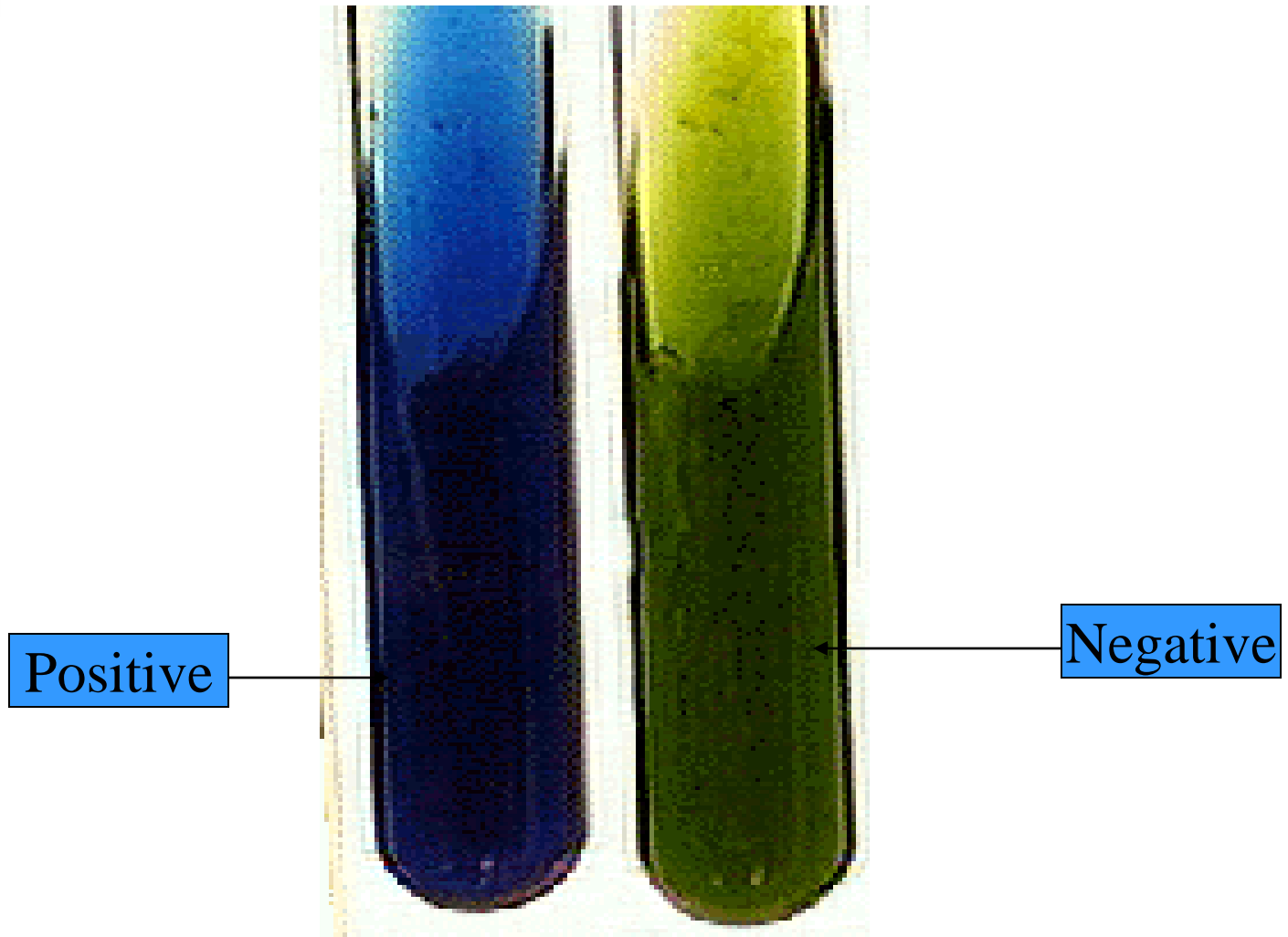

Urea Utilization

- Some bacteria produce **urease**, an enzyme capable of breaking down urea and produce alkaline end products. This distinguishes *Proteus* from other bacteria
- **Media and Reagent:** Urea Broth with phenol red
- **Method:** Inoculate the media with a loop and incubate at 37°C for 24 hours.
- **Expected Results:**
 - **Positive test:** production of alkaline end products = pinkish red color
 - **Negative test:** No color change

Urea Test results

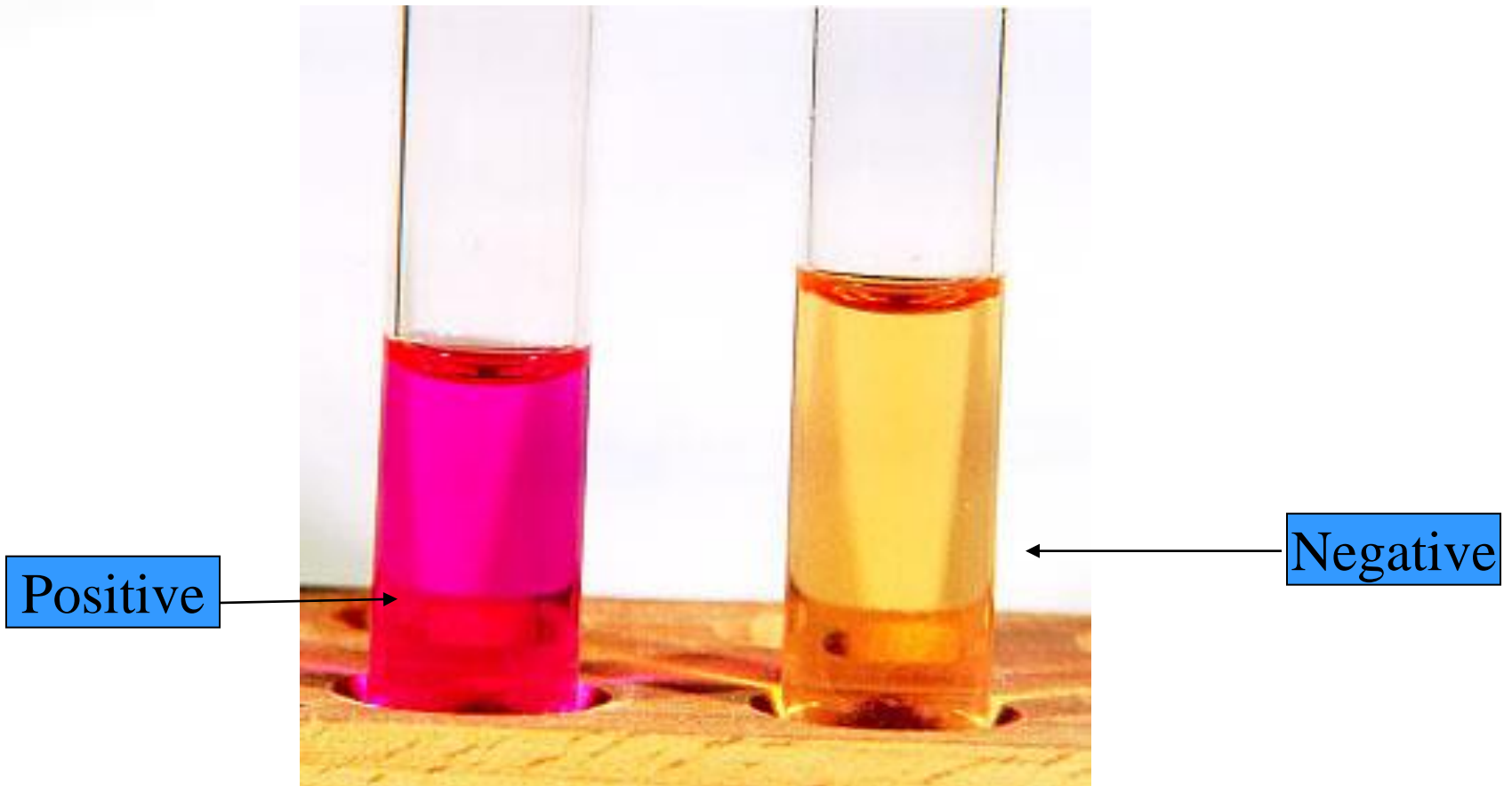


Oxidase Activity

- The oxidase test is a biochemical reaction that assays for the presence of **cytochrome oxidase**.
- In the presence of an organism that contains the cytochrome oxidase enzyme, the reduced colorless reagent becomes an oxidized colored product

Oxidase activity

- The reagent: (N,N,N',N' tetramethyl –p-phenylenediamine) is used to do the test
- The test must be interpreted within 10 to 20 sec, many organisms in this family can give delayed false positive reaction .
- The use of metal loops (due to iron oxide on its surface) to transfer the colonies for the test gives false positive.
Only wooden applicator sticks can be used.

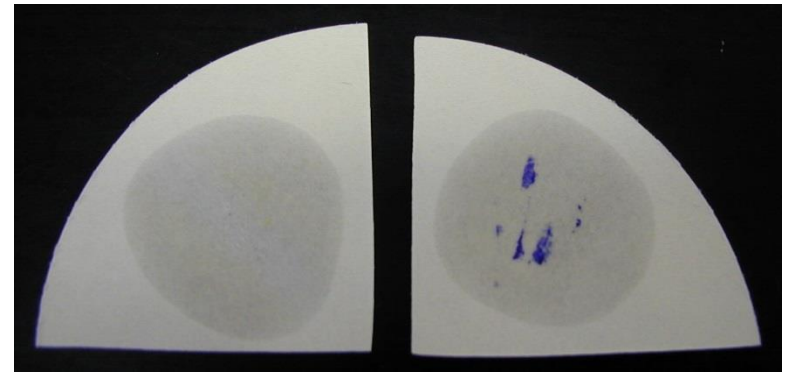

www.mesacc.edu

Nitrate Reduction

This test is used to determine the ability of an organism to reduce nitrate (NO_3) to nitrite (NO_2) using the enzyme **nitrate reductase**.

Nitrate Reduction

- It can be done on any basal medium that support the growth of the organism and contains 0.1% KNO₃.



- Detection of nitrate reduction can be done by:
 - Addition of alpha-naphtyl amine and sulfanilic acid will form diazonium compound (red color).
 - If color did not develop, it means one of two possibilities:
 - NO₃ is not reduced to NO₂ or
 - NO₃ is reduced to NO₂ and further reduced to Nitrogen gas
 - To check for gas production: add a pinch of zinc dust, the development of red color indicates that nitrate reduction did not take place.
- No color development indicates that NO₃ was reduced to NO₂ and further reduced to nitrogen gas

Nitrate Reduction

Nitrate (NO_3)

Nitrate reductase

Nitrite (NO_2)

Nitrite reductase

Nitric Oxide (NO)

Nitric Oxide reductase

Nitrous Oxide (N_2O)

Nitrous Oxide reductase

Nitrogen (N_2)

Negative-Positive Species:


N. mucosa
B. catarrhalis
K. denitrificans



Nitrite (NO_3)

+ Sulfanilic Acid
(Nitrate
Reagent A)

+ alpha-
naphthylamine
(Nitrate
Reagent B)

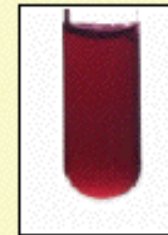

Nitrite (NO_2)

A red color is observed
as Nitrate Reagents A
and B react with
nitrite produced from
nitrate by nitrate
reductase

or

No change in color is
when Nitrate Reagents
A and B are added
because nitrite produced
from nitrate has been
further reduced.



+ Zn dust



Nitric Oxide (NO)
Nitrous Oxide (N_2O)
or
Nitrogen (N_2)

Catalase

- Catalase is an enzyme found in most bacteria. It catalyzes the breakdown of hydrogen peroxide to release free oxygen.
- $2 \text{H}_2\text{O}_2 \text{ -----} \rightarrow 2 \text{H}_2\text{O} + \text{O}_2$
- Procedure: Add one drop of H_2O_2 to a glass slide with a loopful of growth from each culture to be tested. The development of bubbles is indicative of a positive catalase test.
- **The test is performed on a blood-free medium. Why??**



5-5 Catalase test. The test is performed by adding 3% hydrogen peroxide (H_2O_2) to a colony on a glass slide or by adding colony paste on a wooden stick to a drop of H_2O_2 on a slide, as shown here. The appearance of bubbles indicates that the enzyme, catalase, has hydrolyzed H_2O_2 into oxygen plus water. Staphylococci and micrococci are differentiated from other aerobic gram-positive cocci by a positive catalase test (*right*). No bubbles appear in a negative test result (*left*).

Coagulase


- Coagulase is an enzyme that catalyzes the formation of a **fibrin clot in plasma**.
- The presence of coagulase can be detected by heavily inoculating the test organism into rabbit plasma and incubating the mixture for 4 to 24 hours.
- Any degree of clotting during this time, from a loose clot suspended in the plasma to a solid, immovable clot is a **positive** result.
- *Staphylococcus aureus* produces coagulase enzyme while *Staphylococcus epidermidis* does not.

Coagulase Results


- **Reading Results:**
 - If the organism has coagulase it will clump the plasma.
 - If the organism does not have coagulase it will not clump the plasma.



Gram Negative Identification Flow Chart

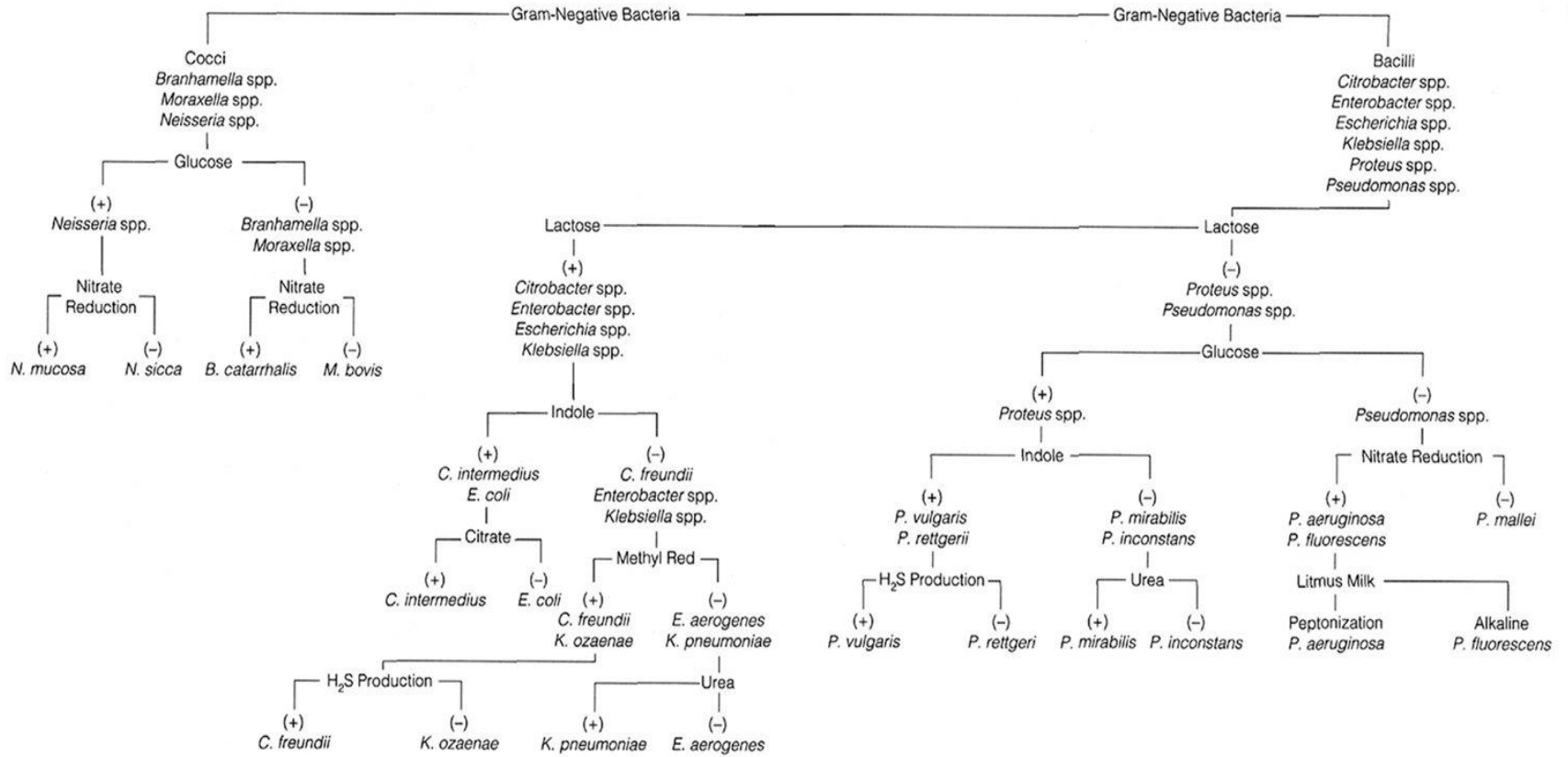

NG: No growth; G: Growth; A/G: Acid and gas; A: Acid only

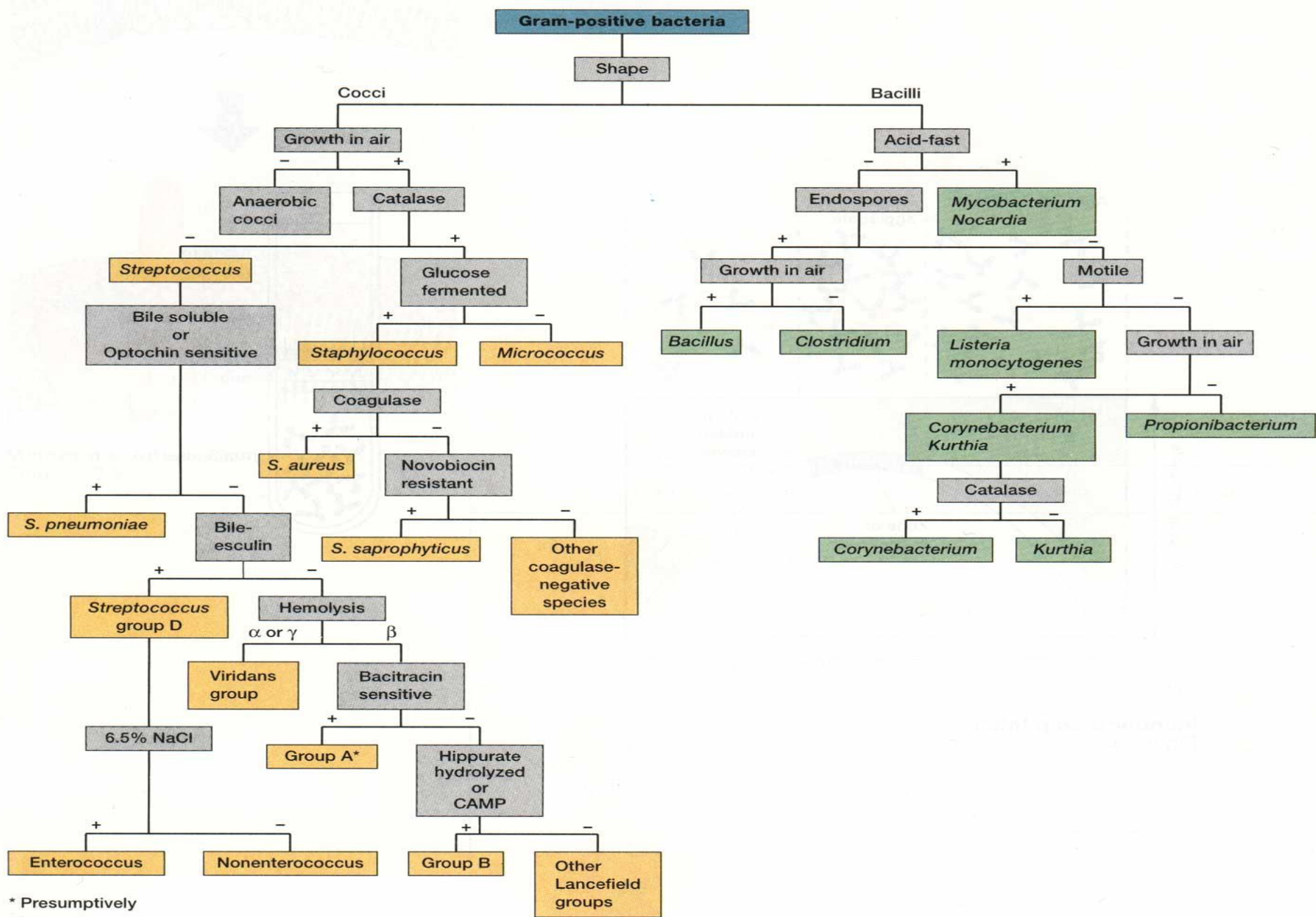

Dichotomous Keys for Clinically Important Genera
 Figure 34.7 a